




MOTORE PER IMPLANTOLOGIA TOUCH SCREEN

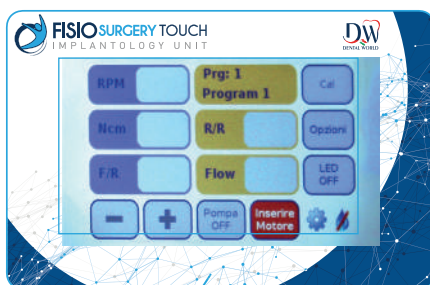

DESCRIZIONE

Fisio Surgery Touch è un'unità per implantologia con schermo touch-screen progettata in modo semplice, sicuro ed affidabile.

La potenza del motore brushless, in associazione alla dualità del torque e della velocità, lo rendono tra le migliori unità presenti sul mercato.

CARATTERISTICHE

- Velocità del motore **da 400 a 40,000 Rpm**
- Valore del torque **da 6 a 72 Ncm**
- 3 differenti flussi di irrigazione: **100 ml / min**
- Rapporto di riduzione:
- **1:1 - 16:1 - 20:1 - 24:1 - 32:1 - 64:1 - 80:1**
- Motore brushless
- Pompa di irrigazione esterna monouso

VANTAGGI

- 5 programmi memorizzabili
- Compensa automaticamente la perdita di torque dovuta al contrangolo

PACKAGING

- Console di controllo elettronica
- Micromotore con luce LED
- Pedale "DeLuxe"
- Asta di sospensione per sacchi da irrigazione
- Tubi per pompe peristaltiche
- Alimentatore



TUBI PER L'IRRIGAZIONE




PEDALE DELUXE

00390097	Fisio Surgery Touch
00390094	Carrying case
00390096	Surgical irrigation tubing - 10 pcs





IMPLANTOLOGY UNIT TOUCH SCREEN

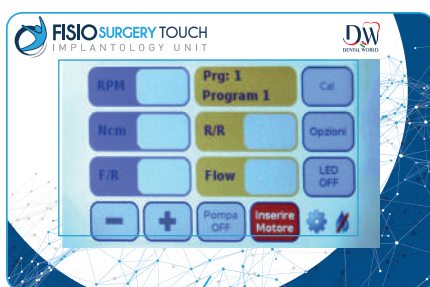

DESCRIPTION

Fisio Surgery Touch is a touch-screen implantology unit designed in a simple, safe and reliable manner.

The power of the brushless motor, in combination with the duality of torque and speed, make it among the best units on the market.

CHARACTERISTICS

- Speed of the motor **from 400 to 40,000 Rpm**
- Torque value **from 6 to 72 Ncm**
- 3 different irrigation flows: **100 ml / min**
- Reduction ratios:
- **1:1 - 16:1 - 20:1 - 24:1 - 32:1 - 64:1 - 80:1**
- Brushless motor
- External disposable irrigation pump

ADVANTAGES

- 5 storable programs
- Automatically compensates for loss of torque due to the contra angle

PACKAGING

- Electronic control console
- Micromotor with LED light
- "DeLuxe" foot pedal
- Suspension pole for irrigation bags
- Tubes for peristaltic pumps
- Power supply



IRRIGATION TUBES



DELUXE FOOT PEDAL

00390097	Fisio Surgery Touch
00390094	Carrying case
00390096	Surgical irrigation tubing - 10 pcs

